
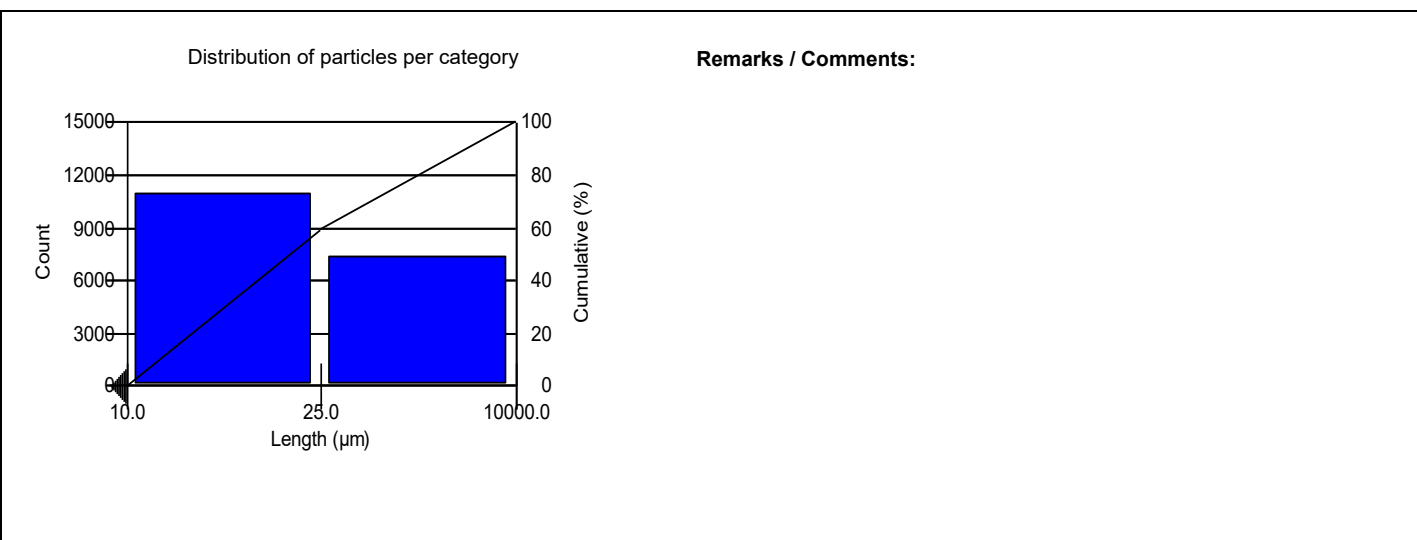


Particulate Matter in Ophthalmic Solutions as per USP 789 (2007)

Date:	2026-May-28, 8:57:36 AM GMT-04:00	Covered area: 1103.15 mm ²	Illumination type: incident
Analyst:	Clemex Support (2026-Ma ...)	Covered ratio: 6.24	Magnification: 55.5x *
Sample ID:		Effective filtration diameter: 15 mm	Calibration: 1.4451 µm/pixel
Quantity of parts:	1.0	Effective filtration area: 176.71 mm ²	Units: Microns
Sample volume:	300 ml		# Fields: 731

Size range (mic)	10 to 25 microns		25 to 50 microns		50 microns and over		
Average number of particles allowed	50	/ml	5	/ml	2	/ml	
Blank count	(can be filled in manually or linked to an exported Vision file)						Total
Particle counts	10965		7310		0		18275
Estimated particles for the filter / piece	1756.5		1171.0		0.0		2927
Corresponding Particle count per ml	5.85		3.90		0.00		9.8
Result of test:	Passed		Passed		Passed		Passed

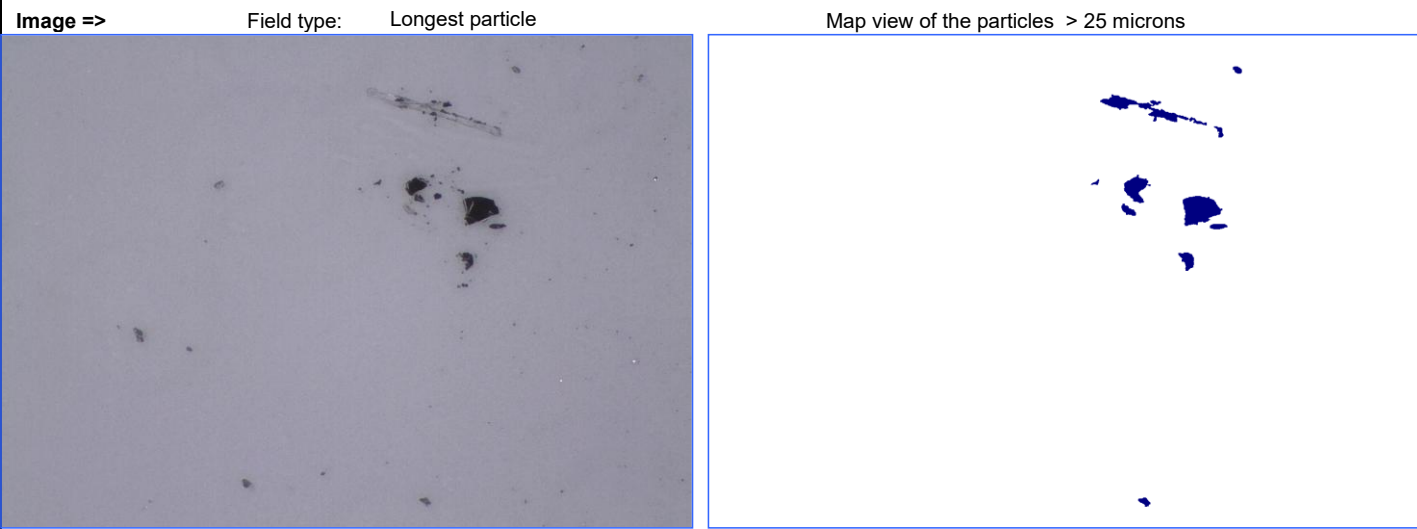

* Note: USP 789 specifies to use 100 +/- 10x magnification.